

GE Healthcare

A grayscale ultrasound image of breast tissue, showing several large, rounded, anechoic (dark) areas, likely representing cysts or fluid-filled structures, surrounded by echogenic (brighter) tissue. The image is used as a background for the top half of the advertisement.

Look differently.

Invenia ABUS

Automated Breast Ultrasound



Invenia™ ABUS

from GE Healthcare offers a
view beyond mammography,
with breast screening technology
that looks differently.



The unseen risk.

Having dense breasts increases
a woman's likelihood to develop cancer

four to six times.¹



Which makes it challenging to detect breast cancer
with mammography alone.

This is a significant problem. Over 40% of women
in the United States have dense breast tissue,² which
can mask the appearance of tumors and limit the
performance of mammography. As breast density
increases, the accuracy of mammograms decrease.

¹Boyd, et al. NEJM Jan 2007.

²Tabár L, et al. Swedish two county trial: impact of mammographic. 2011;260:658-63.



Too many tumors are missed.

For a woman with dense breasts,
screening with mammography
just isn't enough.


Mammography may miss over 1/3 of cancers in dense breasts.³ When it is warranted to look beyond mammography, there now is a screening technology that can effectively detect cancer to deliver confidence and peace of mind.

³Mandelson et al. Breast density as a predictor of mammographic detection: comparison of interval- and screen-detected cancers. J Natl Cancer Inst 2000; 92:1081-1087.



Enhanced detection.

Imaging dense breasts with the
Invenia ABUS can find cancers
mammography misses.

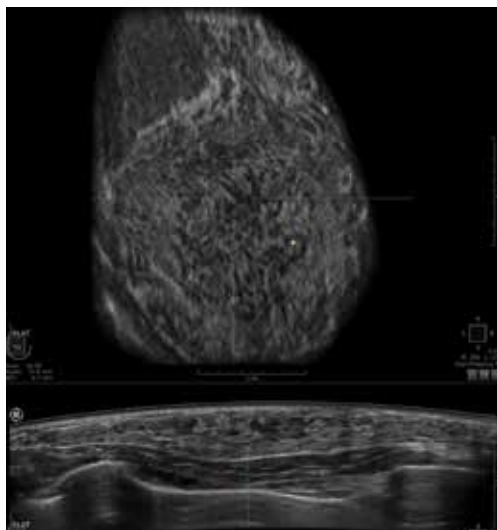

GE Healthcare is the sole provider of FDA-approved technology designed for screening women with dense breast tissue. Invenia ABUS has been shown to improve breast cancer detection by 35.7% over mammography alone.⁴

Multiple clinical research studies⁴ demonstrate that when used as an adjunct to mammography, Invenia ABUS can detect invasive, node-negative cancers.

Lateral view of right breast:

Navigation with the coronal plane efficiently highlights potential abnormalities and streamlines the screening workflow. Superior detail and contrast resolution provides exceptional visualization within dense breast tissue.

⁴Source: FDA PMA P110006 summary of safety and effectiveness.



Designed for screening.

From its remarkable architecture
to its advanced imaging algorithms,
Invenia ABUS is engineered for
automated screening.

Powerful Imaging Architecture

Invenia ABUS's imaging architecture shifts traditional ultrasound from hardware- to software-based processing, resulting in extraordinary performance for the fast-paced breast imaging environment. With its massive parallel processing power and proprietary beamforming technology, the system creates focus at every pixel, delivering an image of high uniformity and resolution.

Intelligent Imaging Algorithms

Advanced algorithms automate the imaging process to help provide remarkable image quality and reproducibility from user to user.

Tissue Equalization Algorithm uses local gradient information to enhance edges and overall system gain, especially within the dense tissue of the nipple and areola complex.

Nipple Shadow Compensation adjusts the signal level under the nipple to be similar to the surrounding area for enhanced visualization.

Breast Border Detection estimates the boundary between breast tissue and the background to mask out non-tissue, high-noise areas – especially important for 3D and coronal imaging.

Chest Wall Detection algorithm detects the chest wall surface based on the breast tissue recognition to eliminate pixels under the surface in 3D. In combination with breast border detection, only pixels in the breast tissue area are used for visualization.

User-friendly touchscreen

Operator workflow is smooth and easy with the Invenia ABUS high-resolution touchscreen display with advanced Projective Capacitive Touch (PCT) technology. Its sleek, graphical user interface enhances the way you work. Tap and swipe colorful icons to quickly and easily maneuver through the Invenia ABUS exam.



Shaped for a woman's anatomy.

Gently **curved** to follow the natural contour of the breast.



The Reverse-Curve™ transducer enhances both patient comfort and breast coverage during the exam.

The 15 cm wide field-of-view transducer automatically creates uniform compression across the entire breast for consistent, reproducible image quality independent of the operator.

The convergent scan line geometry minimizes beam refraction. The result is deep penetration, sharp focus, resolution at depth and enhanced anatomical detail of complex breast tissue and structures.




Compression Assist

Automatically applies compression to the breast for patient comfort, operator ease and image acquisition quality.


A comfortable and comforting exam.

Invenia ABUS helps deliver greater clarity and peace of mind for patients and providers alike — all via a simple, automated procedure that provides a comfortable, non-ionizing alternative to other supplemental screening options.

The evidence is in the images.

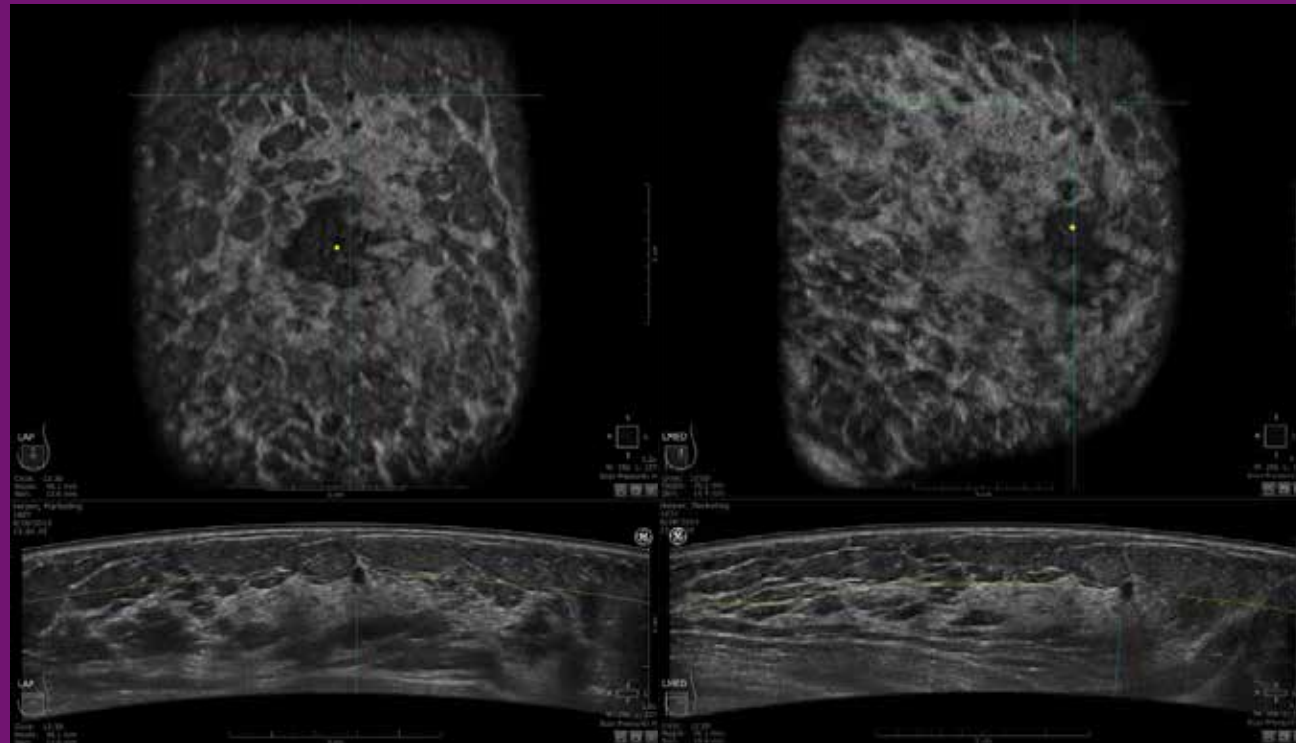

Multi-planar and multi-volume correlation between the left anterior-posterior and left medial views allows for efficient and concise confirmation of this ill-defined mass.

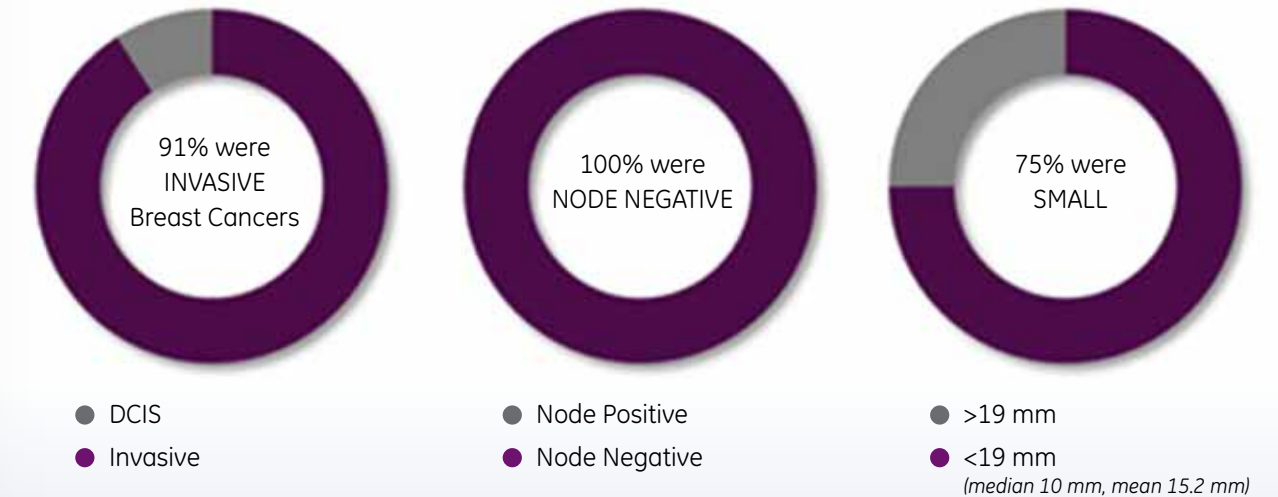

Multiple clinical research studies demonstrate that radiologists can detect more cancers at an earlier, more treatable stage when using ABUS as an adjunctive screening tool with mammography. Ninety-one percent of ABUS-detected cancers with normal or negative mammography were pathologically confirmed as invasive, in this asymptomatic screening population of women with dense breast tissue and no prior breast interventions. The majority of these invasive malignant lesions were also small and node negative.


Study results compiled from USI, 20082002, clinicaltrial.gov NCT00816530 data

A fast, efficient flow.

Designed for fast, efficient breast ultrasound workflow, Invenia ABUS Workstation displays 3D volumes in a patented, 2-mm-thick coronal slice from the skin to the chest wall. You can review, quickly interpret and archive patient exams.

Leveraging the 3-megapixel high resolution monitor, the Invenia ABUS Workstation provides intelligent tools for analysis and for efficient reading and reporting workflow. The system is DICOM 3.0 compliant and integrates with other imaging information systems and repositories.

The Invenia ABUS Workstation generates a customizable, DICOM.pdf format report. Captured images, annotations, and 3D lesion location (clock position, distance from skin and from nipple) are automatically included.

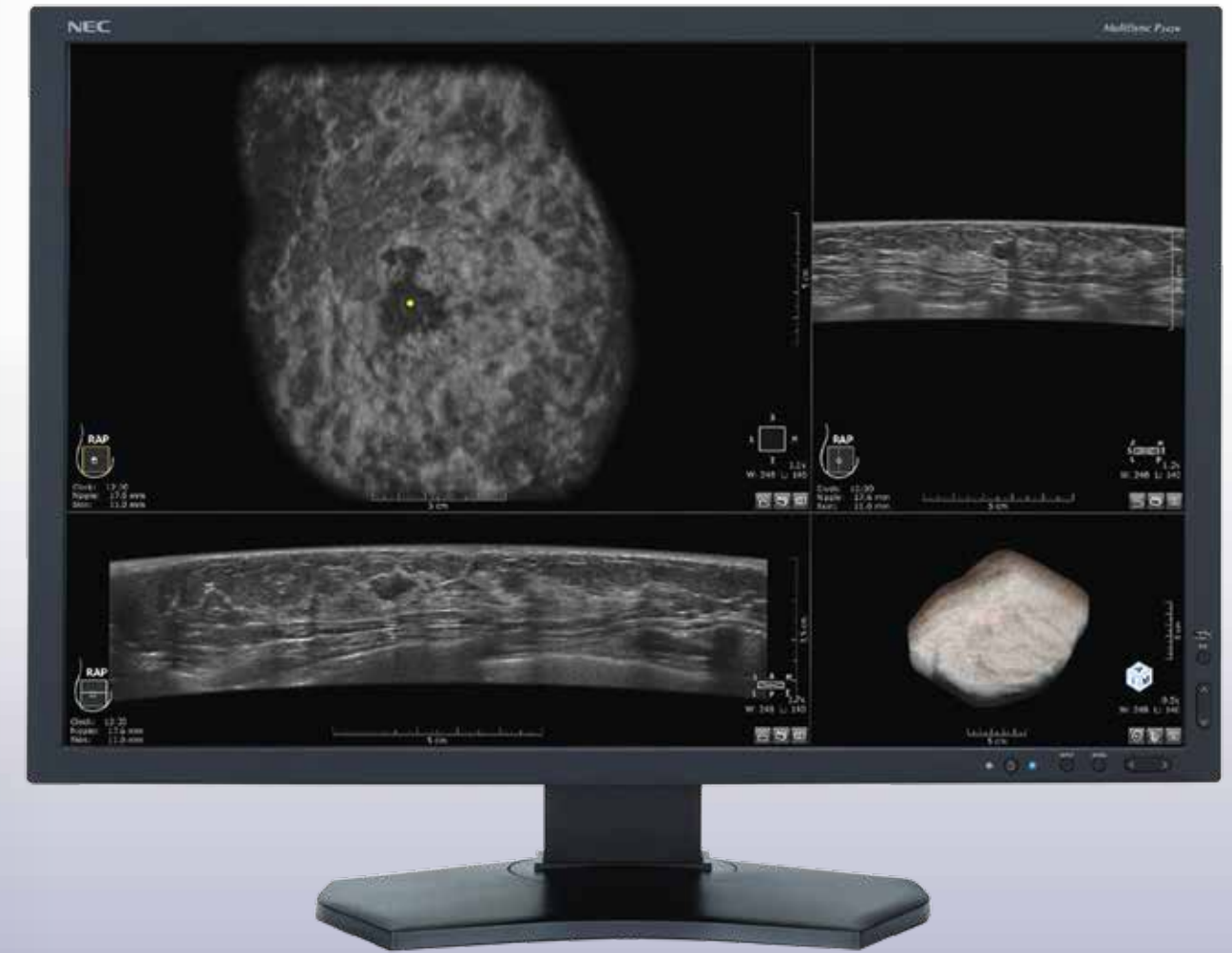



Support beyond technology.

Customizable marketing tools

Targeted marketing tools feature clinically researched and professionally designed content to help you promote Invenia ABUS screening to referring physicians, their patients, and women in your community. Use the customizable templates to inform referring physicians about how Invenia ABUS improves early cancer detection for their patients. Increase women's knowledge about breast density with patient education tools such as brochures, posters and videos. Market to your community with social media, press releases and other templates to help grow awareness.

Comprehensive education and training

The Invenia ABUS Mastery Program uses progressive teaching techniques led by our experienced team of peer educators and clinical applications specialists. The program

is dedicated to the rapid development of diagnostic confidence to quickly build the clinical skills and experience necessary to effectively use Invenia ABUS.

Comprehensive service solutions

Our dedicated field service engineers, online engineers and online applications specialists make it their mission to maximize your asset performance. Additionally, with the InSite™ EXc option, users have access to a suite of service support tools that allow for remote diagnostics and technical assistance.

To support your total cost-of-ownership goals, a variety of flexible service offerings are available.



Growing awareness.

For the 4 out of 10 women with dense breasts, screening with Invenia ABUS provides an additional option for breast care.

Invenia ABUS can find the cancers mammography misses, supporting early detection for better patient outcomes.





Invenia ABUS

Look differently.



About GE Healthcare

GE Healthcare provides transformational medical technologies and services to meet the demand for increased access, enhanced quality and more affordable healthcare around the world. GE (NYSE: GE) works on things that matter – great people and technologies taking on tough challenges. From medical imaging, software & IT, patient monitoring and diagnostics to drug discovery, biopharmaceutical manufacturing technologies and performance improvement solutions, GE Healthcare helps medical professionals deliver great healthcare to their patients.

Better health for more people

Healthymagination is a GE shared commitment to create better health for more people. GE is investing in and innovating technologies that address doctors' and patients' needs. Invenia ABUS was created to improve the quality of breast care.

To learn more about Invenia ABUS, visit gehealthcare.com/ultrasound, call 888-202-5528 or scan the QR code.



healthymagination



imagination at work

Brief Statement

The Invenia ABUS is indicated as an adjunct to mammography or breast cancer screening in asymptomatic women for whom screening mammography findings are normal or benign (BI-RADS® Assessment Category 1 or 2), with dense breast parenchyma (BI-RADS Composition/Density 3 or 4), and have not had previous clinical breast intervention. The device is intended to increase breast cancer detection in the described patient population. The Invenia ABUS may also be used for diagnostic ultrasound imaging of the breast in symptomatic women. See the device manual for detailed information, contraindications, warnings, precautions, potential adverse events.

©2014 General Electric Company — All rights reserved.

GE, the GE Monogram, Invenia, Reverse Curve, and InSite imagination at work are trademarks of General Electric Company.

GE Medical Systems Ultrasound & Primary Care Diagnostics, LLC, a General Electric company, doing business as GE Healthcare.

DICOM is a registered trademark of the National Electrical Manufacturers Association. BI-RADS is a trademark of the American College of Radiology. Third party trademarks are the property of their respective owners.

GE Healthcare
9900 West Innovation Drive
Wauwatosa, WI 53226
USA

April 2014
ULT-0562-01.14-EN-US
DOC1436228